

Tobii Pro Eye Tracker Manager

User's Manual



Version 1.12
01/2019
All rights reserved.
Copyright © Tobii Pro AB (publ)

The information contained in this document is proprietary to Tobii Pro AB. Any reproduction in part or whole without prior written authorization by Tobii Pro AB is prohibited. Products that are referred to in this document may be either trademarks and/or registered trademarks of the respective owners. The publisher and the author make no claim to these trademarks. While every precaution has been taken in the preparation of this document, the publisher and the author assume no responsibility for errors or omissions, or for damages resulting from the use of information contained in this document or from the use of programs and source code that may accompany it. In no event shall the publisher and the author be liable for any loss of profit or any other commercial damage caused or alleged to have been caused directly or indirectly by this document.

Content subject to change without notice.


Please check the web site www.tobii.com for updated versions of this document.

tobiipro

Tobii Pro Eye Tracker Manager is a free configuration and setting utility for Tobii Pro Eye Trackers.

Application Menu
You can find information about the application, such as the current version number and its name.

Connected Eye Trackers
The list of connected Eye Trackers is displayed here, the currently-selected Eye Tracker is highlighted. Disconnected Eye Trackers are grayed out.

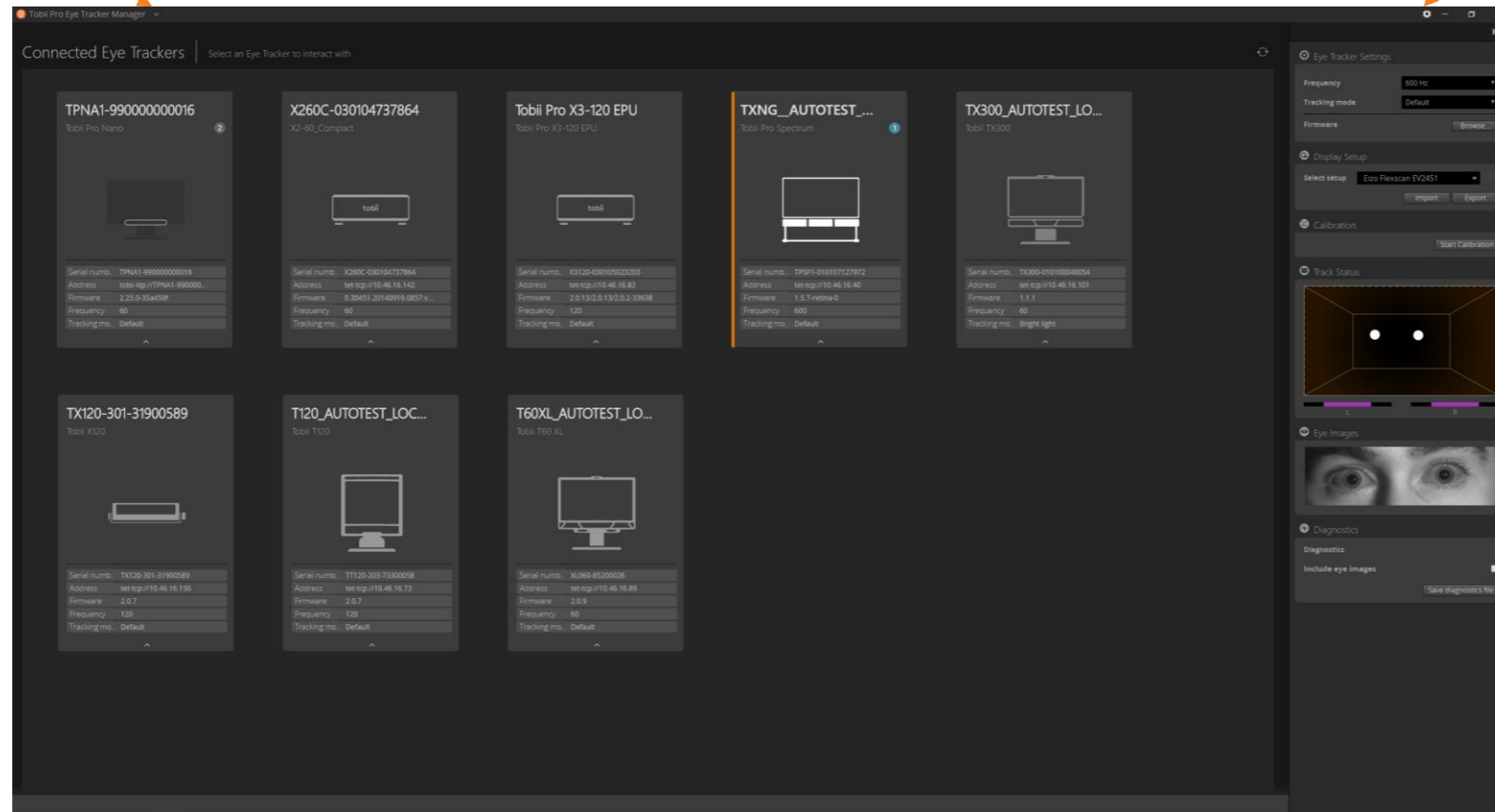

Selected Eye Tracker Panel
You can change settings for the selected eye tracker. Hide or show this panel with the arrow icons (▶▶ or ◀◀).

Application Settings
You can change language and notification settings here.

Eye Tracker Settings
You can change your Eye Tracker's setting such as frequency and tracking mode here. Some of these features are available for specific Eye Trackers only.

Display Setup
You can configure the Eye Tracker with a screen or a physical scene here by creating a new display setup or using an existing one.

Calibration
You can calibrate the Eye Tracker for the user's eyes by clicking start calibration and following the instructions.

Position Guide
The user can find the best position to sit in front of the eye tracker using the Position Guide. Move your head to align with the reference circle for optimal tracking position.

Eye Images
The live view of eye images is shown here. This feature is available only for Tobii Pro Spectrum.

Diagnostics
You can download diagnostic file for support issues here. This feature is available for specific Eye Trackers only.

System Requirements

Operating System	Windows 64 bit & 32 bit	Linux 64 bit	Mac OS X
Versions	7, 8.1 and 10	Ubuntu 16.04	10.11.6 (El Capitan) & 10.12 (Sierra)
CPU	Intel Quad Core i5 3.0 GHz	Intel Celeron 700 MHz	Intel Quad Core i5 3.0 GHz
Memory	2 GB RAM	512 MiB RAM	8 GB 1600 MHz DRAM

Supported Eye Trackers

	Firmware Update	Eye Tracking Mode	Eye Tracking Frequency	Position Guide	Display Setup	Calibration	Diagnostics	Eye Images
Tobii Pro Spectrum	•	•	•	•	•	•	•	•
Tobii Pro TX300		•	•	•	•	•	•	
Tobii Pro T60XL				•		•	•	
Tobii X60		•		•	•	•	•	
Tobii T60		•	•	•		•	•	
Tobii X120		•	•	•	•	•	•	
Tobii T120		•	•	•		•	•	
Tobii Pro Nano				•	•	•		
Tobii Pro X3-120				•	•	•	•	
Tobii Pro X2-60				•	•	•	•	
Tobii Pro X2-30				•	•	•	•	